

Automatic Load and Unload System for AMSCO™ 600 Steam Sterilizers



Taking Productivity
to **NEW HEIGHTS**

REDEFINING WORKFLOW

RECLAIM PRECIOUS TIME with the Automatic Load and Unload System (ALUS) for the AMSCO™ 600 Steam Sterilizer



LIFTING PRODUCTIVITY TO NEW HEIGHTS

Focus on other time-sensitive requests while instruments are automatically loaded, processed and unloaded.

EASY CHAMBER ACCESS

Hinged side shields and a removable table make it easier to access test packs and the chamber for cleaning or repairs.



KEEP YOUR PROCESS IN MOTION

FLEXIBILITY & SAFETY



Configure for auto load/unload
or manual load/auto unload

Hinged side shield = extra safety

Table sensor = proper load position

Automatic or manual
operation modes


Casters lock table in place

EFFICIENCY

- Reduce stress and strain of transporting heavy loads and decrease time spent loading the sterilizer.
- Chamber pass-through allows return of empty car to loading area.

AUTOMATIC CYCLE SELECTION

Queue up your next load using the bar code reader to select your cycle in advance.

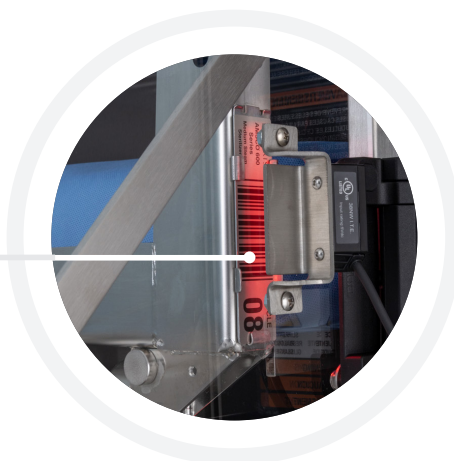

IT JUST FITS

The sleek design of the **AMSCO 600 Medium Steam Sterilizer** meets your needs – maximizing productivity in a minimum footprint.



The Automatic Load and Unload System (ALUS) for the AMSCO™ 600 Steam Sterilizer

Available Sizes	26.5" x 26.5" x 39" Length (673 mm x 673 mm x 991 mm)	26.5" x 26.5" x 51" Length (673 mm x 673 mm x 1295 mm)	26.5" x 26.5" x 63" Length (673 mm x 673 mm x 1600 mm)			
Available Configurations	Door Options	# Table	Manual Load	Auto Load	Auto Unload	Barcode Reader
	Double Door	1	✓		✓	
	Double Door	2		✓	✓	
	Double Door	2		✓	✓	✓
	Single Door	1		✓	✓	

Contact your STERIS representative for more information.