

Central Sterile Services

DISCOVER
POWERFUL
WASHING

AMSCO™ 7000HP SERIES
WASHER-DISINFECTOR

**PERFORMANCE,
PRODUCTIVITY
and PRECISION**



One Integrated Approach to Healthcare


 **STERIS**



PERFORMANCE

Optimal cleaning with complete rack coverage


- Ensure full coverage with cross-pattern spray arm with 120° range
- High and low impingement options

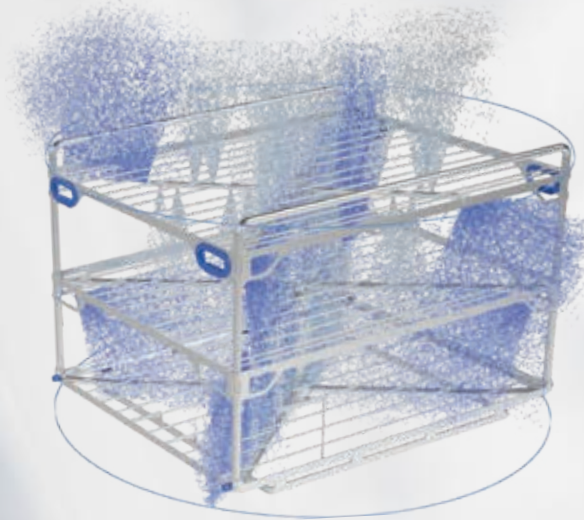

Dislodge stubborn soils

from cannulated devices


The 3-, 4- and 5-Level Manifold Racks can be equipped with up to 10 lumen flushing hoses

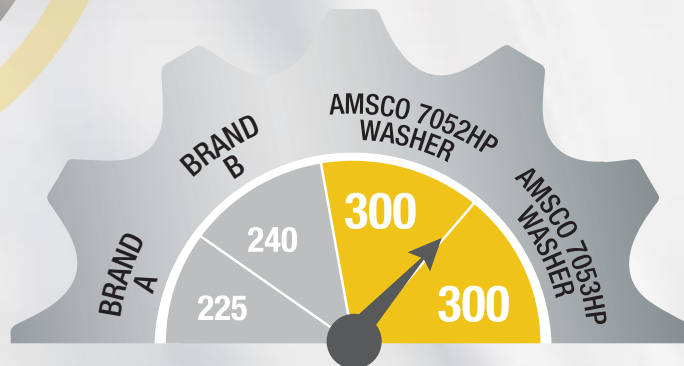

DIN Trays per 8 hour day

If trays don't fit correctly, how can spray arms clean them properly?

With other competitors' racks, it doesn't matter how you load DIN trays...

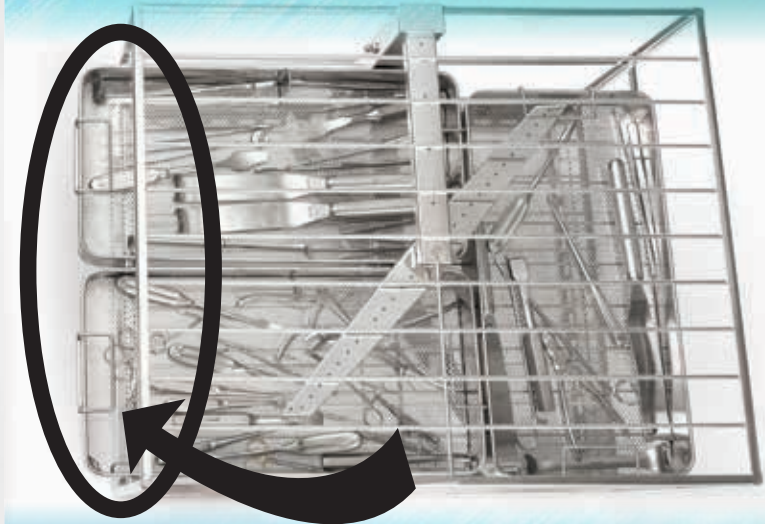

...they don't fit, but with **STERIS** racks they do

Unlock the **industry's fastest cycle time** with the power of Prolystica™ High Performance (HP) Cleaning Chemistries



Deliver maximum cleaning efficacy

with precise chemistry dosing from integrated flow meters

Increase productivity and reduce washer-disinfector maintenance cost with low-foaming formulations:

- Maintain pump pressure, reducing possible cycle aborts and related service costs
- Maximise spray arm performance and mechanical soil removal

Powerful combination

STERIS considered **each and every factor** in your process, from the instruments to the equipment processing them, to develop formulations that **positively impact your cleaning processes**

Prolystica HP Cleaning Chemistries provide superior instrument and capital protection

Automated Advantage

Prolystica HP Cleaning Chemistries reduce unnecessary wear and tear in your washer-disinfector – **providing savings beyond the cost of chemistry**





PRODUCTIVITY



AMSCO 7000 Series HP Washer-Disinfectors take productivity to elevated heights, offering **the industry's highest throughput***

Stay on schedule and informed

with cycle progress notification lighting



Washer-disinfector available



Cycle in progress



Requires attention

Easy cycle selection

with icon-based control panel



The simplicity of connectivity

Seamlessly connect and get cycle information from your AMSCO 7000 HP Washer-Disinfectors with ConnectAssure Technology

Whether data is sent to your facility's tracking system or to ConnectAssure's interface, your cycle information is available with the click of a mouse for essential compliance reporting




*When compared to current healthcare washer-disinfector models



Continuous throughput

using cool-touch handles to remove loads immediately after cycle completion

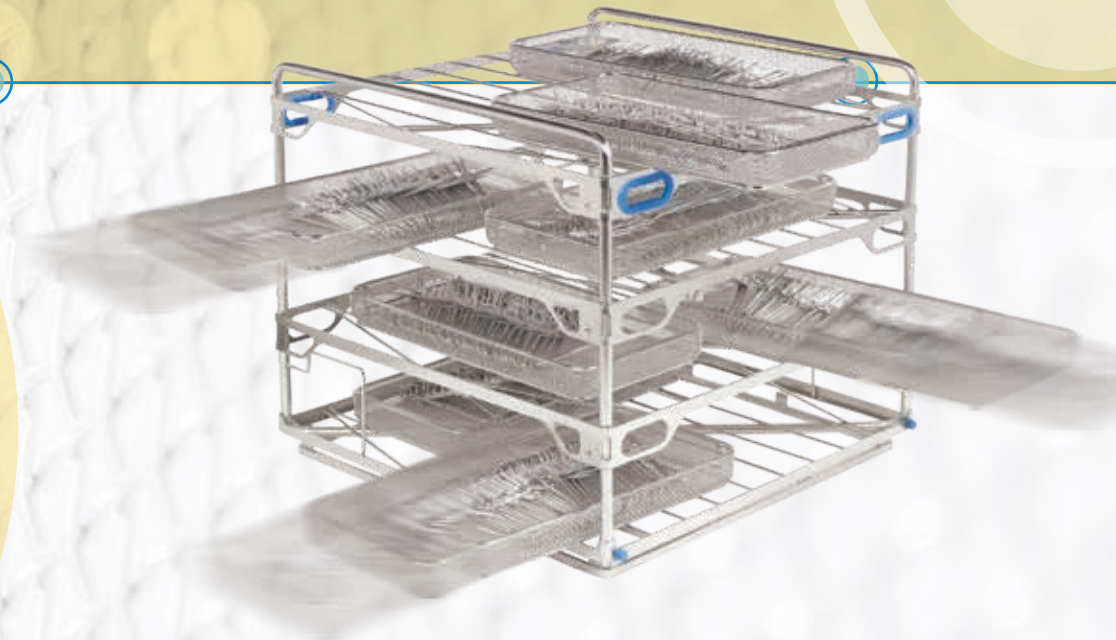

Quicker load / unload preparation

with the versatility of 360° loading, combined with cool-touch handles for immediate safe handling

Flexibility to process a wide range of surgical instruments with

multiple rack configurations





SUSTAINABILITY WITH PRECISION

Combining fast, effective cycle times with reduced water consumption minimises environmental impact while reducing utilities usage



AMSCO 7052HP Washer-Disinfector uses up to

50% LESS

water per cycle*

Acu-Fill

Reduce water usage by **4** litres

and improve start up times by **30** seconds

Acu-Wash

Recover and reuse **11** litres

of thermal rinse water in every cycle

*Sources: Performance Specifications relative to STERIS and competitive brands

UP TO **23** litres & NEEDED **3** minutes

Acu-Rinse

Reduce water consumption by

8 litres

and cycle time by

2 minutes

Acu-Dry

Reduce cycle time by


30 seconds

or more when temperature sensors detect and complete cycles as soon as dry

AMSCO 7053HP Washer-Disinfecter uses up to

**40 %
LESS**

water per cycle*

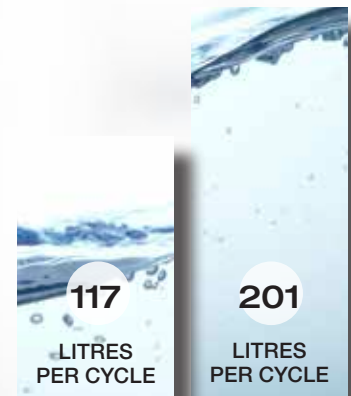

AMSCO
7053HP

Brand A

TOTAL SAVINGS PER CYCLE!

7052HP Washer-Disinfector

Overall external dimensions (WxHxD):	1067 x 2051 x 813 mm (142 x 80-3/4 x 32")
---	--

7053HP Washer-Disinfector

Overall external dimensions (WxHxD):	1067 x 2051 x 965 mm (142 x 80-3/4 x 38")
---	--



EMEA CLINICAL EDUCATION
STERIS
UNIVERSITY

Professional education & clinical expertise

Get connected to STERIS University, encompassing training disciplines and quality clinical education, delivered by a team of expert industry professionals.

For more details, ask educationdesk@steris.com or visit university.steris.com

 **STERIS**

STERIS

Chancery House
Rayns Way
Leicester LE7 1PF
United Kingdom
Tel: +44 116 274 0600
Email: contact_hc@steris.com
www.steris.com

